


**voco**

# HIGHLIGHTS

**High-Performance-Composite**  
GrandioSO Unlimited –  
4 mm universal composite  
in five Cluster-Shades

**Vital and suitable for daily use**  
Preservation of pulp vitality and  
pulp capping with the biomaterial  
MTA vpt (Vital Pulp Therapy)

**Tasteful**  
Profluorid Varnish+BioMin  
with bioglass is the new power  
for enamel hardening

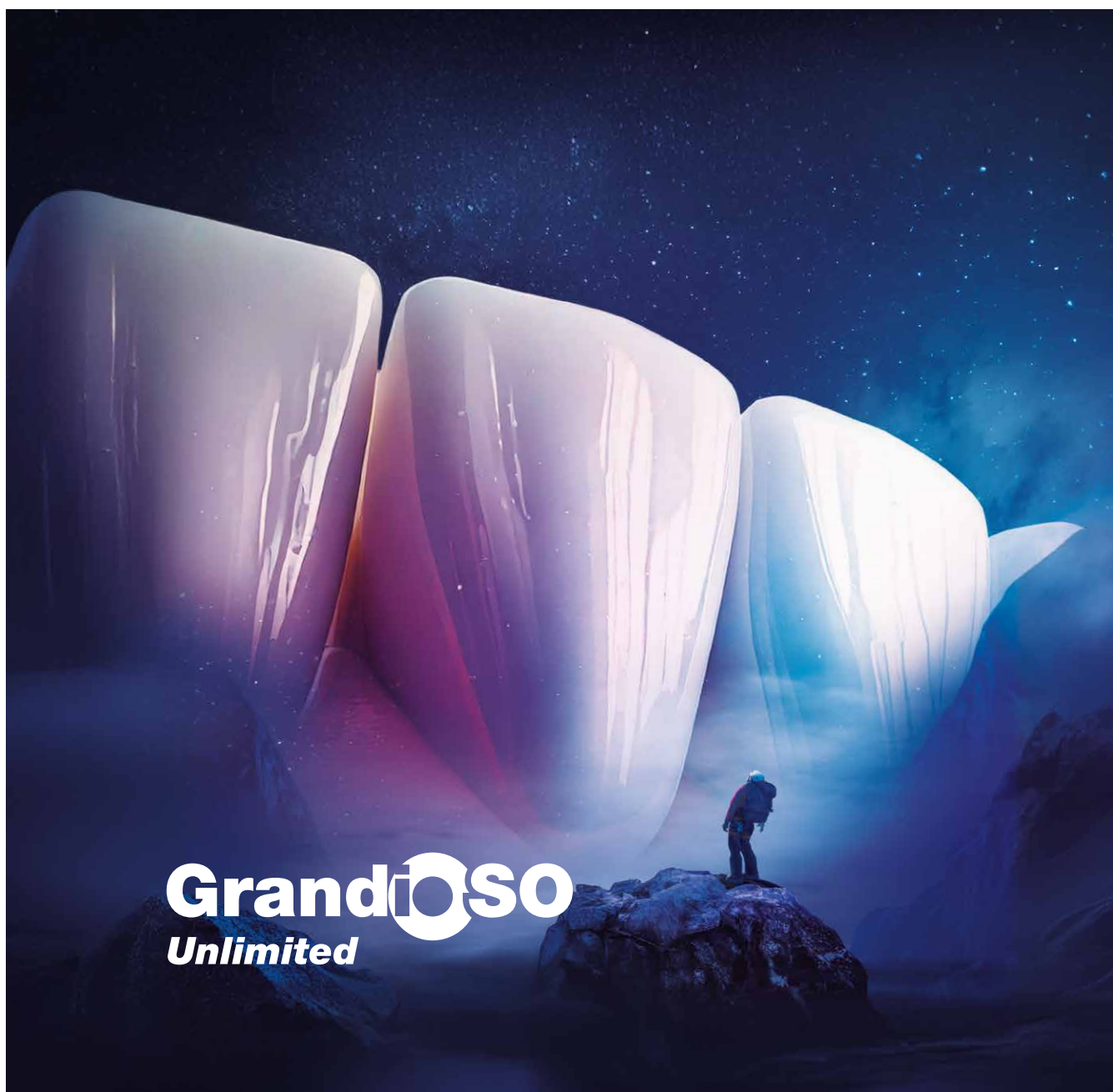

**GrandioSO**  
*Unlimited*

# GrandioSO Unlimited

4 mm universal composite in five Cluster-Shades

While aesthetic universal composites have to be layered and polymerised in 2 mm increments, which is complicated, bulk-fill composites offer the advantage of application into the cavity in increments of up to 4 mm. However, due to their insufficient aesthetics, bulk-fill composites are normally only used for the posterior region. VOCO has now resolved this issue with the new universal composite GrandioSO Unlimited. The aesthetic material can be applied in increments of up to 4 mm into the cavity and is indicated both for the anterior and posterior regions. You therefore only need one packable universal composite!

## Cluster-Shade-System

In addition, with only 5 shades to stock, you are optimally positioned for your everyday practice and can significantly downsize your stock. Thanks to the simplified shade system, the 5 Cluster-Shades are able to cover all 16 classic VITA® shades:

GrandioSO Unlimited	A1				A2		A3			A3.5				A4		
VITA® classical	A1	B1	C1	D2	A2	B2	A3	B3	D4	A3.5	C2	C3	D3	B4	A4	C4

## High-Performance-Composite

With an impressive filler content of 91 % by weight, GrandioSO Unlimited is setting new standards in the high-performance composites class when it comes to stability and strength. The low volume shrinkage of 1.44 % reduces the danger of poor marginal gap formation and marginal discolouration, and contributes to the most high-quality restorations.



## Test pack

10 × 0.3 g Caps  
 (2 each of × A1, A2, A3, A3.5, A4)  
 + 10 × Futurabond U SingleDose,  
 + shade card

REF 6257



## Syringe set

5 × 4.5 g syringe  
 (A1, A2, A3, A3.5, A4)  
 + shade card

REF 6258

## Caps set

75 × 0.3 g Caps  
 (15 each of × A1, A2, A3, A3.5, A4)  
 + shade card

REF 6270

Further presentations and shades can be found under [www.voco.dental](http://www.voco.dental)

## The High-Performance-Composite

**Universal**  
For anterior and posterior regions

---

**Highest stability and strength**  
91 % w/w filler content

---

**5 Cluster-Shades**  
Covering all 16 VITA® shades

---

**No layering**  
Up to 4 mm increments

---

**Excellent aesthetics**  
Outstanding polishability

---


# MTA vpt

Maintaining pulp vitality in daily practice –  
biocompatible and monomer-free

A biocompatible, stable and bacteria-tight pulp seal is essential to preserve the vitality of the pulp. With MTA vpt, VOCO offers a Mineral Trioxide Aggregate (MTA) material that has been specially designed for this purpose. Within the scope of Vital Pulp Therapy (VPT), MTA vpt can be used for cavity lining, indirect and direct pulp capping as well as for the treatment of the pulp following vital amputation (partial or complete pulpotomy).

As part of caries lesion management, residual carious dentine can also be reliably sealed with MTA vpt.

MTA vpt releases calcium and hydroxide ions during the setting phase, which stimulate the body's own regeneration processes and the formation of tertiary dentine. The resulting high pH is known for its bacteriostatic effect and supports remineralisation. In this way, MTA vpt creates a stable and bacteria-tight protective barrier – the optimum basis for maintaining the vitality of the pulp.

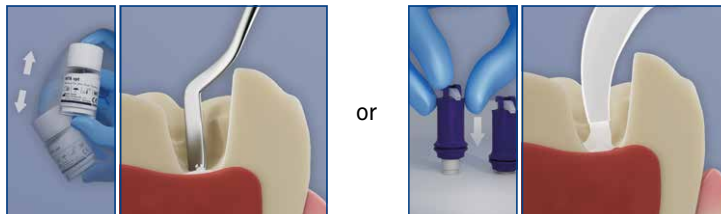

- **Reliable:** mineralises to form a stable bacteria-proof seal
- **Time-efficient:** immediate further processing thanks to layering without pressure
- **Practical:** for use in selective or stepwise caries management and vital pulp therapy



**MTA vpt**  
handmix  
10 g powder + 10 ml liquid  
REF 1645



**MTA vpt AC**  
25 pcs. Application capsule  
REF 1646

**Avoid Endo.**  
**Save the pulp.**

Preserving pulp vitality in daily practice

releases  
OH<sup>-</sup> OH<sup>-</sup>  
Ca<sup>2+</sup>

MTA

Vital Pulp Therapy

290% radiopacity



Fluorescence that misses nothing:

## Bifix Temp fluorescent

Dual-curing, composite-based temporary luting material

5 ml QuickMix syringe  
+ accessories

REF 1056



Fluorescent luting material

Handpiece, battery charger,  
power cord, 2 device bases,  
rubber shield, protective sleeves,  
protective glasses

REF 9070



## AlignerFlow LC

Fluorescent composite for aligner attachments

2 × 2 g syringe A1

REF 1526

2 × 2 g syringe A2

REF 1527

2 × 2 g syringe A3

REF 1528




Fluorescent composite for aligner attachments

CURE & CHECK

## TwinWave

- High performance: State-of-the-art with high output of 1,300 mW
- Lightweight and handy design
- Easy to use: Mode (light polymerization) can be selected using a practical button
- Continuously ready for use: Rechargeable battery

detects residues, saves time, ensures safety

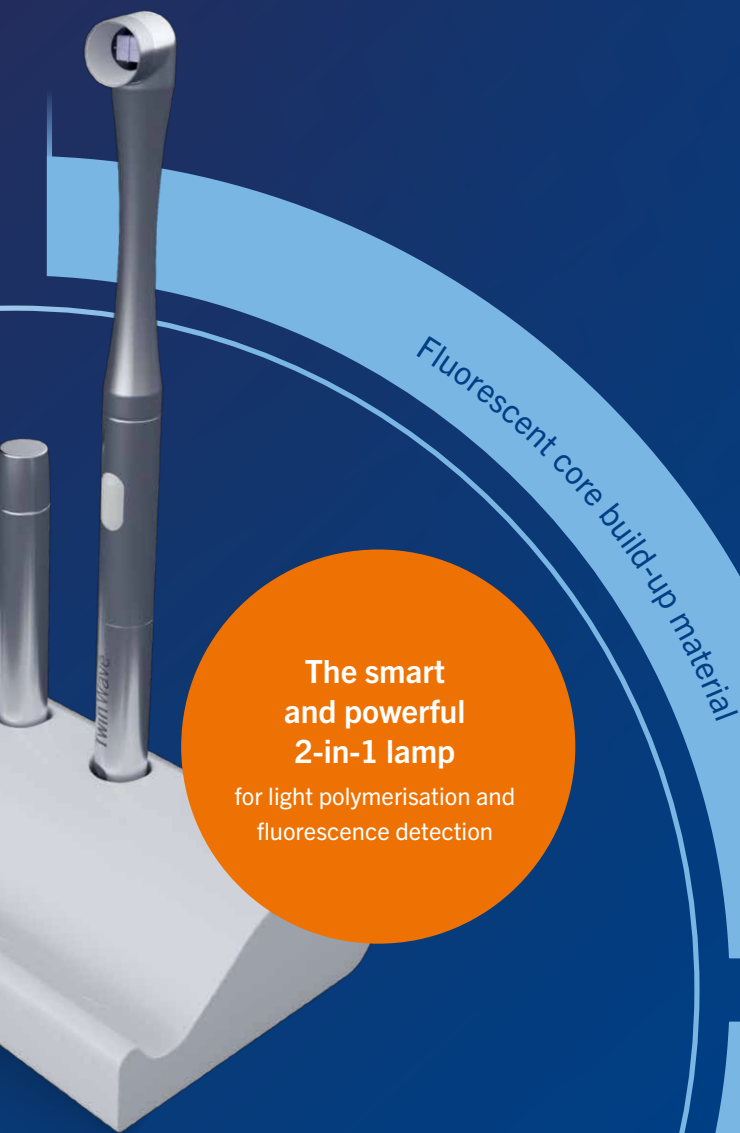

The smart  
and powerful  
2-in-1 lamp

for light polymerisation and  
fluorescence detection

ave

**NEW**

ne-art LED technology  
/cm<sup>2</sup>

: Metal housing, total weight only 75 g

nerisation, fluorescence detection)

cal rotary switch

replacement battery included

Fluorescent filling material for temporaries

## Rebilda DC fluorescent

The first dual-curing flowable core build-up and  
post luting system with activatable fluorescence

50 g Cartridge,  
fluorescent dentine,  
+ accessories

REF 1399

10 g QuickMix syringe  
fluorescent dentine,  
+ accessories

REF 1406



## Clip F fluorescent

Light-curing filling material  
for temporary restorations

2 × 4 g syringe

REF 1288



# Profluorid Varnish + BioMin

Save  
**25%\***

Significantly hardens weakened enamel –  
in four pleasant flavours

The established fluoride varnish – VOCO Profluorid Varnish – gets the addition of bioglass from BioMin. This invisible bioglass dissolves in the presence of saliva and releases calcium and phosphate ions, which support the remineralisation of the tooth. In combination with the fluoride from the varnish, crystalline structures form on the tooth surface, which significantly harden acid-weakened enamel, for example.

A real double protection together with the well-known desensitisation!

## 200 × 0.40 ml SingleDose

Caramel REF 2276  
Cherry REF 2277  
Mint REF 2278  
Bubble Gum REF 2279

## 50 × 0.40 ml SingleDose

Caramel REF 2270  
Cherry REF 2271  
Mint REF 2272  
Bubble Gum REF 2273



## 10 ml Tube

Caramel REF 2266  
Cherry REF 2267  
Mint REF 2268  
Bubble Gum REF 2269

## Mixed Flavour

48 × 0.40 ml SingleDose  
sorted: caramel, cherry,  
mint, bubble gum  
REF 2274

# Profluorid Varnish

Fluoride-containing dental desensitising varnish (5 % NaF)

Hypersensitivity is mainly caused by exposed dentinal tubules. Therefore, the aim of treatment is to seal these tubules. Profluorid Varnish is a fluoride-containing varnish to effectively seal dentinal tubules exposed on the necks of teeth or after cavity preparation.

## The advantages at a glance

- Fast desensitisation and fluoride release (22.600 ppm fluoride)
- Easy to apply – moisture tolerant
- Available in seven pleasant flavours
- Economical application in thin layers, high yield



Save  
**25%\***

## 5 × 1,70 ml Cartridge

Melon REF 2235

## 10 ml Tube

Melon REF 1267 Caramel REF 2232  
Mint REF 2234 Cherry REF 2233  
Pina Colada REF 2237 Cola Lime REF 2236

## 200 × 0.40 ml SingleDose

Melon REF 1271 Bubble Gum REF 2239  
Mint REF 2231 Cherry REF 2230  
Pina Colada REF 2243 Cola Lime REF 2241  
Caramel REF 2229

## Mixed Flavour

4 × 10 ml tube  
sorted: mint, melon,  
cherry and caramel  
REF 2971

# Futurabond U

In the practical *SingleDose*

Behind Futurabond U lies an extensive research and development process that had one goal in mind: to offer a universal adhesive, giving the practitioner nothing but advantages and reliability and thus eliminating the need for any other adhesive systems in the dental practice.

Futurabond U sets new standards in regard to bond strength, and especially with its flexible and versatile application.

## 50 *SingleDose*

50 × Futurabond U  
+ accessories

REF 1571

## 200 *SingleDose*

200 × Futurabond U  
+ accessories

REF 1572

Save  
**12%\***



# Grandio blocs & Grandio disc

Nano-ceramic hybrid CAD/CAM-material  
for permanent restorations and veneering of frameworks

Grandio blocs/Grandio disc impress with their physical values in terms of elasticity, strength and a high filler content of 86 % by weight. As a result, it is closer to the natural tooth than any other material for indirect restorations. Grandio blocs/disc is the material of choice and is superior to conventional restoration materials. The products can be easily finished intraorally and polished to a high gloss.



Available in  
12 mm, 15 mm  
and 18 mm

## Grandio disc multicolor

1 × 15 mm (ø 98 mm)

For example:

A2 REF 6217

A3 REF 6218



Source: Dr. Felipe Moura, Brazil

Further presentations and shades  
can be found under [www.voco.dental](http://www.voco.dental)

Also  
available for  
Planmill®



## Grandio blocs

5 × no. 12 LT/HT

REF see table

## Grandio blocs

5 × no. 14L LT/HT

REF see table

	For the anterior region							For the posterior region			
	A1 LT	A2 LT	A3 LT	A3.5 LT	B1 LT	C2 LT	BL LT	A1 HT	A2 HT	A3 HT	A3.5 HT
No. 12	6003	6004	6005	6006	6007	6008	6009	6012	6013	6014	6015
No. 14L	6018	6019	6020	6021	6022	6023	6024	6027	6028	6029	6030




# Bifix Veneer LC

Safe and aesthetic –  
veneer luting with optimised workflow

The particularly aesthetic light-curing nanohybrid luting composite Bifix Veneer LC filled with 81 % by weight was specially developed for the simple and safe adhesive luting of veneers, inlays and onlays made of composite or ceramic.

The clear shade concept with the **three shades bright, clear and warm** makes placing veneers easy and enables a positive effect on the aesthetic result. Thanks to the precisely coordinated try-in pastes, choosing the right shade is easy and the final result can be checked in advance.

**An intuitive and efficient shade concept for maximum aesthetic success**

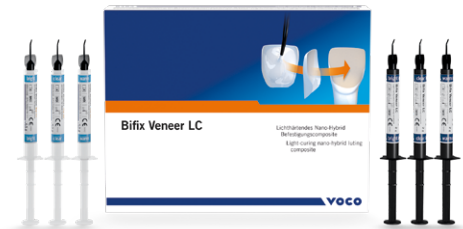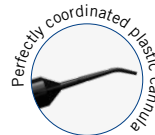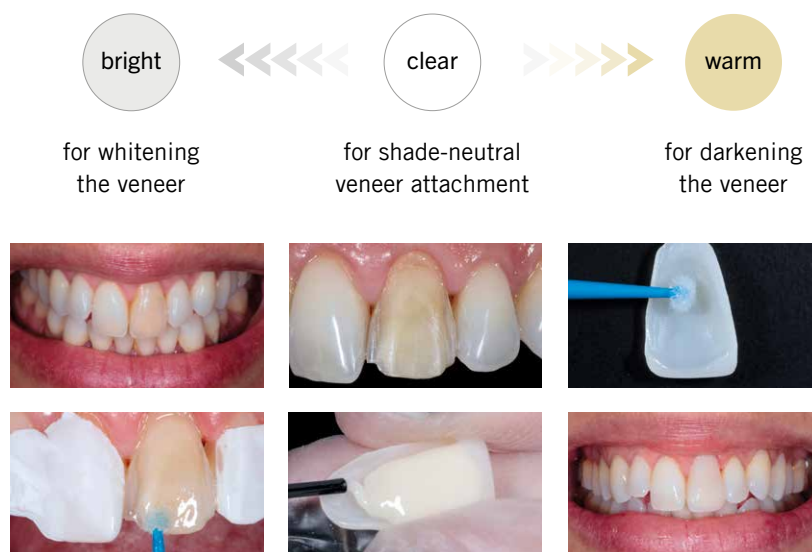

## Bifix Veneer set

3 × 2 g syringe Bifix Veneer LC  
3 × 2.5 g syringe Bifix Veneer Try-In  
(bright, clear, warm)  
+ accessories

REF 1428



## Bifix Veneer System set

3 × 2 g syringe Bifix Veneer LC  
3 × 2.5 g syringe Bifix Veneer Try-In  
(bright, clear, warm)  
2 ml syringe Vococid  
20 pcs. *SingleDose* Futurabond U  
20 pcs. *SingleDose* Ceramic Bond  
10 pcs. Sticky Tim  
+ accessories

REF 1429

Further presentations and shades  
can be found under [www.voco.dental](http://www.voco.dental)

Please contact your dental depot or your VOCO dental consultant  
to order the VOCO Highlights presented here.

\*compared with individual purchase

Offers are only valid at participating retail partners. Special prices offered by your depot are always taken into consideration. Errors and omissions excepted. No liability for printing errors.

**VOCO GmbH**  
Anton-Flettner-Straße 1-3  
27472 Cuxhaven  
Germany

